



26 CLINICAL EVALUATORS

976 TOTAL USES

96% CLINICAL RATING

Category: Disposables | Occlusion Articulating Paper

Description

TrollFoil is a super-thin, double-sided articulating foil designed to provide precise and consistent occlusal markings on wet, dry, and highly polished surfaces. The 8-micron foil is pre-mounted in a semi-rigid frame that keeps it flat and taut, allowing convenient one-handed placement without forceps while preventing ink from transferring to the clinician's fingers. A removable protective strip supports clean handling, and the thin design is suitable for challenging clinical situations, including limited opening, smooth implant surfaces, gag-sensitive patients, and metal or nonmetal restorations.

Indications

- Checking occlusion on wet or dry teeth
- Use in patients with limited mouth opening
- Marking occlusion on implants with smooth surfaces
- Use with gag-sensitive patients
- Checking contacts on metal and nonmetal restorations

Unique Attributes

- Built-in tab allows placement without articulating forceps.
- Foil material coated with carbon-based color, offering more accurate results.
- Ultra-thin foil provides precise, clearly defined contact markings.
- Produces visible marks on Zirconia, Porcelain, Metal, Composite and other restorative materials.
- Performs on both wet and dry surfaces, with red and blue options for distinguishing occlusal movements.



Photo courtesy of Dr. Sabiha Bunek



Clinical Tips

- Dry the tooth surfaces before use to improve the clarity of the markings.
- Use blue foil for centric occlusion and red foil for excursive movements to help identify interferences.
- Have the patient bite firmly, then gently pull the foil through the contact to enhance the markings.
- Make sure the patient's teeth are biting only on the marking paper and not on the frame holding it together.

"THIS IS A
10 OF 10.
THE BEST
WE'VE TRIED."

Evaluators' Comments

"Very thin, so it marks only where the contact points are. It does not tear easily."

"It is a necessary product in all offices."

"The plastic form makes handling it much easier than loose strips."

"The markings are very distinct and easy to read. Point contacts are the clearest I've ever seen with any brand of articulating medium."


"For occlusal equilibrations and adjustments of dental orthotics, I liked it a lot and will continue to use it."

Consultants who would
want to stock in office:

58% Yes, in addition to
current product

35% Yes, instead of
current product

Evaluation Summary:

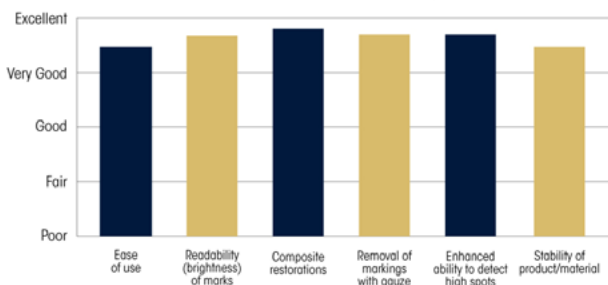

Compared to Competitive Products:

