


THE DENTAL ADVISOR™

Improving Patient Care Through Research & Education

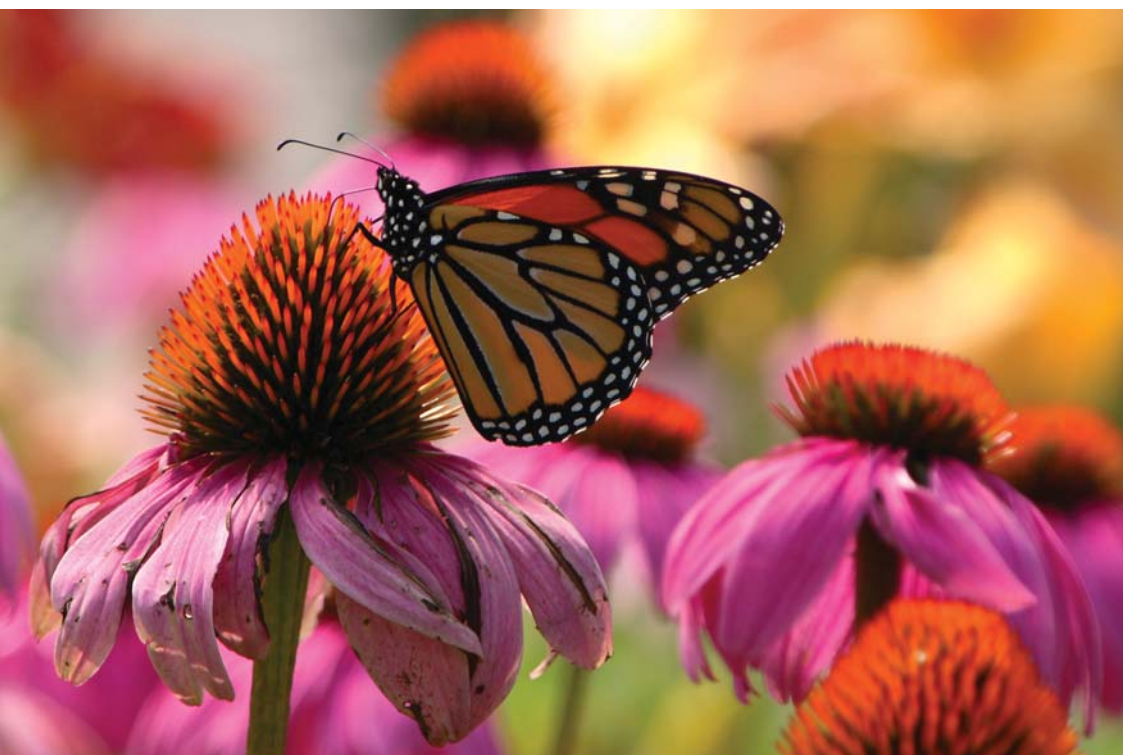

Preventive Hygiene

A dental professional's ultimate challenge is to treat patients in a preventive manner regardless of systemic conditions, compliance with home care, or compromised dentition. This issue of THE DENTAL ADVISOR will focus on specific products found to be helpful by clinical teams in preventive hygiene.

May 2015

Vol. 32, No. 04

MAIN TOPIC

Preventive Hygiene 2

CLINICAL EVALUATIONS

NTI-tss Plus, NTI-tss Plus Soft
(Occlusion - TMJ Device)..... 7

Versa-Temp 2
(Provisional Crown and Bridge Material).... 10

AirBlade AirMirror 10

Crystal HD® Mouth Mirrors 12

LONG-TERM EVALUATION

ACTIVA BioACTIVE-RESTORATIVE
6-month Clinical Performance..... 8

EDITORS' CHOICE

Nupro White
Varnish 6

Vera by Young
(Disposable Prophyl Angle) 11



www.dentaladvisor.com

RATINGS:

Excellent + + + + +

Very Good + + + +

Good + + +

SENIOR EDITORS

John W. Farah, D.D.S., Ph.D.
John M. Powers, Ph.D.

EDITOR-IN-CHIEF

Sabiha S. Bunek, D.D.S.

CLINICAL EDITOR

Lori K. Brown, D.D.S.

EDITORIAL BOARD

John A. Molinari, Ph.D.
Peter Yaman, D.D.S., M.S.
William A. Gregory, D.D.S., M.S.
Santine Harlock, D.D.S.
Julius E. Bunek, D.D.S., M.S.
Brent Kolb, D.D.S.
Nizar Mansour, D.D.S., M.S.
Charles I. McLaren, D.D.S., M.S.
Kathy O'Keefe, D.D.S., M.S.
William T. Stevenson, D.D.S.
Robert J. Stevenson, D.D.S.
Gary Bloomfield, D.D.S.
Eric Brust, D.D.S., M.S.
Michelle Elford, D.D.S.
Stacy Griffith, D.D.S.
Doug Hock, D.D.S., M.S.
Mark Zahn, D.D.S. M.S.

CONTRIBUTING AUTHOR

Ashley Hood, B.S., R.D.H.

EXECUTIVE TEAM

Matt G. Cowen, B.S.
Jackie Farah, M.A.Ed.
Heidi L. Graber
Tricia G. Hurtt
Jennifer Kalasz
Callie Knapp, B.A.
John A. Molinari, Ph.D.
Peri D. Nelson, B.S.
Courtney Richardson, B.F.A.
Christopher Voigtman
Nelson Williams, M.S.
Mary E. Yakas, B.A., CMC
Ron Yapp, M.S.

PUBLISHER

Dental Consultants, Inc.

Please send inquiries and address changes to:
THE DENTAL ADVISOR,
3110 West Liberty, Ann Arbor, MI 48103
Call: 800.347.1330 - 734.665.2020
Fax: 734.665.1648
Email: info@dentaladvisor.com
Website: www.dentaladvisor.com

No unauthorized duplication or reprints may be made. Inquiries concerning duplication may be directed to the publisher. Copyright ©2015, Dental Consultants, Inc. All rights reserved. Printed in the U.S.A. (ISSN 0748-4666) by Print-Tech, Inc.

This publication is printed on paper that is 50% recycled and has 25% post-consumer content.


From the Desk of Dr. Bunek, Editor-in-Chief



It's ironic how the stages in the oral cavity follow the circle of life. We start out edentulous and many times return to this phase, or at least partially. Small children regularly present with an increase in caries due to inadequate oral hygiene, often times resulting from a lack of dexterity or dietary issues (sugar consumption). Elderly individuals often face similar challenges; we may see an increase in caries possibly due to a lack of dexterity and many times coupled with xerostomia, secondary to the use of medications. Whatever phase of life a patient presents to us in, there are specific associated needs and challenges.

It is our job as dental professionals to address these needs not only at the time of treatment, but also to arm our patients with the information and tools to take with them into their homes. In this issue, our guest author, Ashley Hood, R.D.H., addresses specific challenges faced by our patients and identifies recommended products.

As always, I welcome your comments and suggestions; you can reach me at drbunek@dentaladvisor.com. Thanks for your continued support and reading!

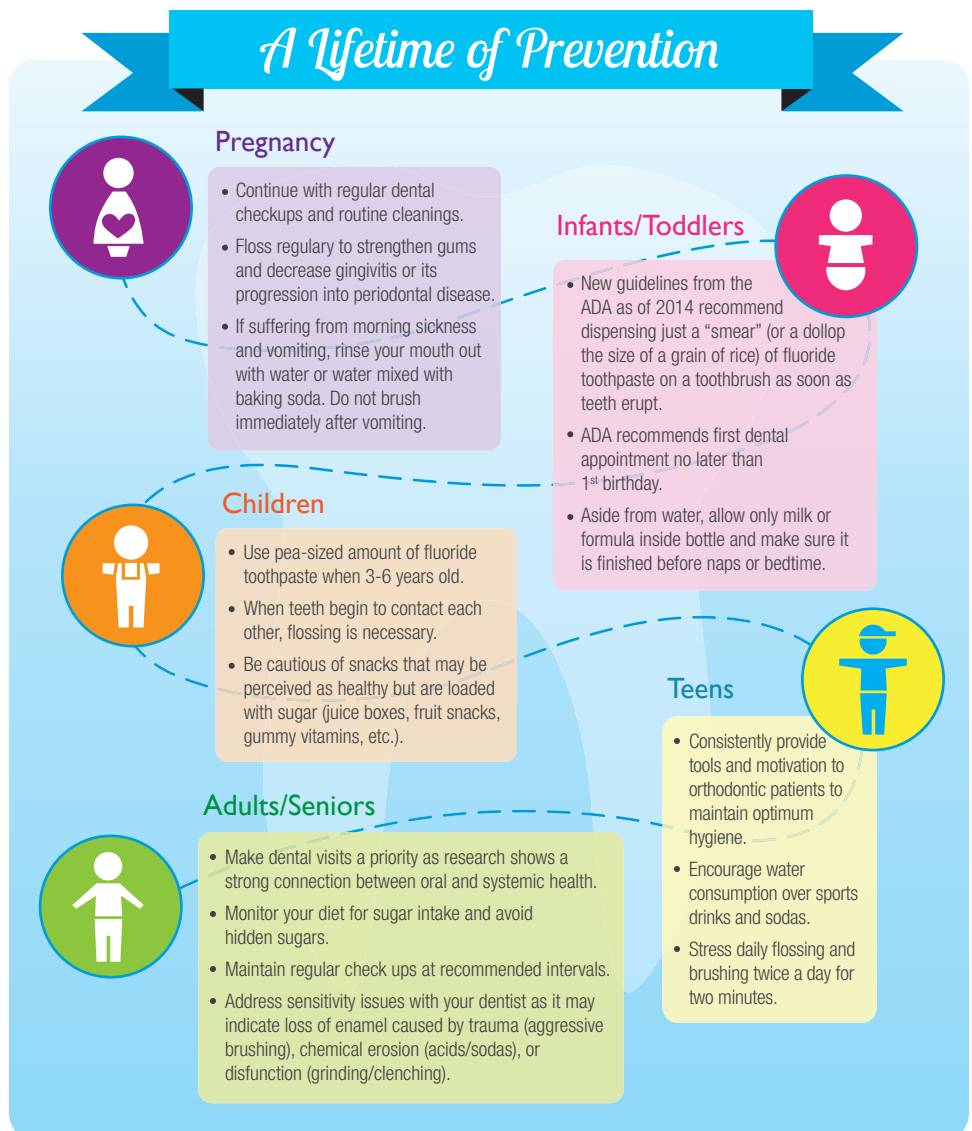

Infants - Teens Dentition Challenges

Challenge: Motivating children to brush and floss



Baby Buddy Oral Care Program
Compac Industries, Inc.
www.smile4babybuddy.com

Baby Buddy is an innovative Oral Care Program, designed to care for the oral cavity as teeth develop. Each stage of 6 stages has an age specific oral care device designed to assist parents and children with a healthy mouth.

Stage 1 (0-16 months) - **Tooth Tissues:** All natural, bamboo wipes with Xylitol. Fluoride free and bubble gum flavored

Stage 2 (3+ months) - **Silicone Finger Brush:** Dual sided, BPA Free silicone, with brush and bumps to massage gums and clean teeth gently as they appear.

Stage 3 (4-16 months) - **Wipe-N-Brush:** a combination of Stage 1 and Stage 2 products, designed for developing teeth.

Stage 4 (4+ months) - **Baby's First Toothbrush:** Designed to be gripped by the infant, 100% BPA free silicone, has soft bristles shaped like an H to clean all teeth at once, while also acting as a teething device

Stages 5 & 6 **360 degree Toothbrush:** The 360 Toothbrush's Bristles are Short, but made of Soft, High Elastic, Ultra Fine Micro Fiber. They are designed to clean all tooth surfaces and reach every corner of the mouth, teeth, gums, and tongue. This allows children to clean using any direction. The stage 5 brush is slightly smaller to accommodate children up to 24 months, and the stage 6 brush is designed to accommodate age 2-12.

Hygienist Comments

"Cloths are very soft for wiping baby teeth."

"Great complete system to teach young parents how to care for the mouth in all stages of development."

"360 toothbrush bristles tend to fall apart fairly easily."

Challenge: Teaching children how to floss properly



Gumchucks Oral Wise
www.gumchucks.com

An innovative flossing concept designed for those with limited dexterity to floss. **Gumchucks** resemble miniature nunchucks, and feature disposable tips that are equipped with a 3/4 inch piece of dental floss. **Gumchucks** are also available for adults (**Gumchucks Adult Pro Floss**) and orthodontic patients (**ORTHOgami**). All products offer an easy way to floss using the C shape recommended by dental professionals.

Hygienist Comments

"Motivates children to floss."

"Allows parent to direct proper flossing technique."

"Innovative app, coloring pages, and characters to keep kids interested."

"Great for oral care education."

"Kids need direction to not be aggressive with flossers as it can damage tissue."

"Costly compared to traditional floss; proprietary handle."

Challenge: Motivating orthodontic patients to maintain optimum oral health



Platypus Orthodontic Flossers
www.platypusco.com

Specially designed orthodontic flossers, ideal for children and adults in traditional orthodontic brackets and wires. Y-shaped flosser allows floss to access difficult areas, and the bracket brush on the opposing end allows patients to clean around brackets.

Hygienist Comments

"Ends the frustration of trying to use traditional floss."

"Keeps mouth clean and free of debris often found in and around wires and brackets."

Innovation in Preventive Care: BrushyBall

Occasionally, we hear of a product that is just a great idea! **BrushyBall** was invented by a dad trying to motivate his daughter to brush. See more information at Brushyball.com

Brushing COMPLETELY is now a follow-the-leader game! **BrushyBall** is designed to motivate young brushers, captivate short attention spans and teach complete brushing in the form of a fun follow-the-leader game. **BrushyBall** speaks for two minutes and music plays while pearly white teeth light up sequentially in six sections to guide children through complete brushing of all teeth and tooth surfaces each time they brush their teeth.



Healthy Adults Dentition Challenges

Challenge: Sensitive Teeth



Remin Pro VOCO

www.vocoamerica.com

Remin Pro is a water based cream containing hydroxyapatite, fluoride, and xylitol indicated for treatment and prevention of hypersensitivity, demineralization and acid erosion. **Remin Pro** is easy to apply and comes in strawberry, melon or mint flavor.

Hygienist Comments

"I love this on my patients post whitening for sensitivity."

"Sugar free."

"Neutral pH."

"Can be used as a tissue protectant and lubricant under dentures."

Challenge: Cleaning around implant crowns



GUM® Soft Picks Sunstar Americas

www.gumbrand.com

Soft Picks are flexible and tapered synthetic rubber bristles that fit easily between most teeth, removing dental plaque where gum disease often starts and toothbrushes alone cannot reach.

Hygienist Comments

"Excellent for cleaning around implants, will not scratch surface."

"Disposable and portable soft picks."

"Cost effective."

Challenge: Providing tools to aid in gum tissue stimulation and removal of debris



Rubber Tip Stimulators Sunstar Americas

www.gumbrand.com

The **GUM® Stimulator** is designed to massage gum tissue and stimulate blood flow. The tips can also be useful for dislodging food particles and clean between teeth.

Hygienist Comments

"Excellent oral hygiene aid for those who find it difficult to floss."

"Excellent for use around orthodontic brackets."

"Tips can be discarded after 3 months and replaced using the same handle."

"Cost effective."

Preventive Success Stories



Submitted by:

Meliza Trevino, RDH
McAllen Family & Sedation Dentistry
McAllen, TX

I would like to share this simple technique I use often to help my patients with generalized or localized recession. I have seen the use of fluoride varnish reduce the sensitivity to my patients more often than not. The feedback I have received is always positive and the patient result speaks volumes. Although varnish is not a covered procedure at times, my patients will opt to pay out of pocket due to the positive result. I follow three steps to placing the varnish; Step one, place Cotton rolls in all areas of treatment, Step two, gently air dry the areas, and Step three, apply the brush on varnish. There is no trick to placing the varnish. I believe localizing the area of sensitivity is key in experiencing a successful outcome.



Submitted by:

Rebecca Dominczak, RDH
Bartoszewicz Family Dental
Grand Rapids, MI

We had a patient that we educated about the link between gum disease and heart disease and recommended daily flossing. The patient went from having several areas of 5 mm pockets and heavy bleeding to having all pocket numbers to 1 to 3mm. The change was incredible!

Senior Population Dentition Challenges

Challenge: Maintaining bone and tissue of endentulous/partially endentulous patients



Nimbus Microfine Brush Nimbus Dental

www.nimbusdental.com

Extra soft toothbrush designed by a periodontist to be gentle on gum tissues. The Microfine brush has two stages of bristles, long tapered to reach difficult areas, and shorter to remove plaque.

Hygienist Comments

"Great for post surgical patients; allows cleaning without further trauma."

"Can use for aggressive brushers who may be causing recession or abrasion. Allows the patient to clean effectively without further damage."

Challenge: Providing an effective alcohol-free mouth rinse to patients

CloSYS Antiseptic Oral Rinse Rowpar Pharmaceuticals

www.closys.com

An alcohol-free, chlorine dioxide rinse for use as a pre-rinse for hygiene appointments, as a post rinse or home care aid for soft tissues.

Hygienist Comments

"Does not stain like chlorhexidine."

"Alcohol free."

Challenge: Increased medications often cause dry mouth

GC Dry Mouth Gel GC America, Inc.

www.gcamerica.com

*Provided in a gel form, **GC Dry Mouth Gel** is an easy delivery for elderly patients or patients with dry mouth due to cancer treatment. The gel can be applied to tissues to minimize pain associated with dry mouth*

Hygienist Comments

"Available in many flavors."

"Sugar free."

"Neutral pH."

"Can be used as a tissue protectant and lubricant under dentures."

THE DENTAL ADVISOR'S Speaker's Bureau



Lecture Topics

The Speakers' Bureau offers continuing education courses on a wide variety of topics for dentists, hygienists, assistants, lab technicians, and office teams. Each topic can be 1-2 hour, half-day and full-day lectures. Topics can be combined and adapted for a custom presentation.

To book a Dental Advisor speaker, or for more information, contact:

Mary Yakas

(734) 665-2020 ext. 129

mary@dentaladvisor.com

Nupro White Varnish

DENTSPLY Professional
(800) 989-8826 | www.dentsply.com

Description

Nupro White Varnish is a 5% sodium fluoride varnish with rapid fluoride release. It is formulated to reduce the amount of time patients must avoid hot beverages, crunchy foods, brushing, etc., by half. **Nupro White Varnish** comes in a unit-dose pack that the user can snap onto a prophylaxis grip, place on a flat surface, or grasp in between fingers. **Nupro White Varnish** can be applied to wet teeth and sets to a clear coating. It is available in 0.25 g and 0.4 g in four flavors: mint (only 0.4 g), bubblegum, raspberry and grape. Each box includes brushes and a prophylaxis grip. **Nupro White Varnish** was evaluated by 36 consultants in 817 uses. This fluoride varnish received a 96% clinical rating.

Suggested Retail Cost	
Box of 50	\$103.00
Box of 100	\$180.00
Box of 500	\$772.00

Product Features

Nupro White Varnish quickly set to a clear finish, which is important to patient acceptance. The variety of flavors was appealing to adults and children. The quick fluoride release that allows patients to be free of restrictions in two hours was well received. Application is smooth, with no clumping or dripping. It can be applied to wet teeth but went on more evenly if the teeth were relatively dry. The **Nupro White Varnish** packs are easy to hold and handle.

Clinical Tips


- Stir varnish well before applying.
- Keep the package flat during application to avoid dripping. ■

+++++




Consultants' Comments
"Patients appreciated the reduced time for restrictions."
"Goes on smoothly."
"Dries clear."
"I would prefer brushes attached to fluoride packet."
"Best varnish I have ever used."


Key Features:

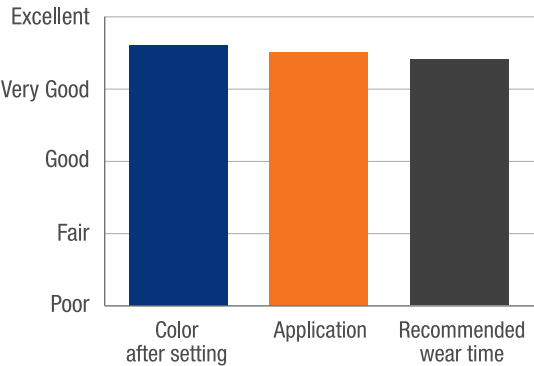

Compared to Competitive Products:

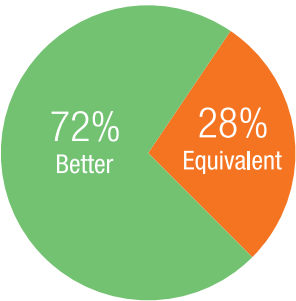

Percentage of Consultants Who Would:

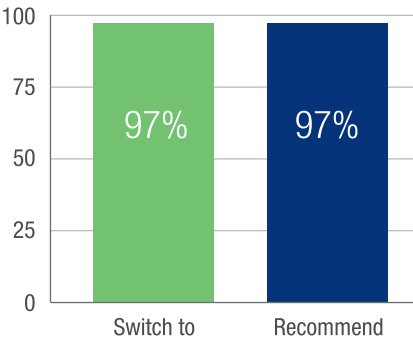

NTI-tss Plus, NTI-tss Plus Soft

Keller Laboratories, Inc.

www.kellerlab.com

Description

NTI-tss Plus is a small acrylic appliance that is worn over the mandibular or maxillary anterior teeth; a discluding element over the incisors prevents occlusion in the posterior. **NTI-tss Plus** is designed to reduce clenching intensity thereby protecting teeth, muscles, and joints while minimizing chair time. It is indicated for therapeutic treatment of the following: tooth wear resulting from bruxism and clenching, migraines and tension headaches, muscle pain associated with parafunction, TMD, and for aiding in diagnostic treatment planning. There are several design options to choose from in order to increase patient compliance. The **NTI-tss Plus Soft** extends coverage to the 2nd bicusps for added retention and is fabricated with a softer inner layer for increased patient comfort. Its indications are the same as above and more specifically for patients with extreme crowding, veneers, worn restorations or natural dentition, and sensitive teeth. **NTI-tss Plus** and **NTI-tss Plus Soft** are laboratory-fabricated appliances. Records from the dentist include accurate maxillary and mandibular models and a measurement of the patient's maximum protrusion. **NTI-tss Plus** and **NTI-tss Plus Soft** were evaluated by 16 consultants who delivered over 150 appliances to patients. These occlusal appliances received an 89% clinical rating.

Product Features

NTI-tss Plus and **NTI-tss Plus Soft** fit accurately and have a nice finish. Delivery appointments are typically very quick since occlusion is limited to the incisors. The appliances were reported to be used for a variety of indications and were most useful for treatment of acute TMD symptoms. **NTI-tss Plus** and **NTI-tss Plus Soft** were reported to be effective in reducing clenching; some patients with chronic headaches reported improvement. The biggest concern of consultants was the potential for unwanted tooth movement, such as anterior open bite, with prolonged use. The discluding element naturally protruded over the lip; patients varied in their tolerance of this. Overall, this small appliance was well accepted by patients.

++ ++



Consultants' Comments

"The lab fabricated them so well that they required little to no adjustment upon delivery."

"The appliance is very effective at relieving my patients' headaches and facial pain symptoms."

*"I feel that **NTI-tss Plus** is great to add to the armamentarium to treat TMJ symptoms."*

"Useful for a patient who is a severe gagger and cannot tolerate coverage of the posterior teeth being."


"Very effective appliance. Patients noted a reduction of symptoms almost immediately."

"When the patients bring their NTIs to their regular appointments, check for changes in occlusion."


Clinical Tips

- If the patient has a severe class II occlusion, make a maxillary appliance.
- Include the maximum protrusive measurement on the RX. Upon delivery, mark the limit of protrusion and trim to furthest mark. This is now the perfect appliance for the patient. ■


Key Features:

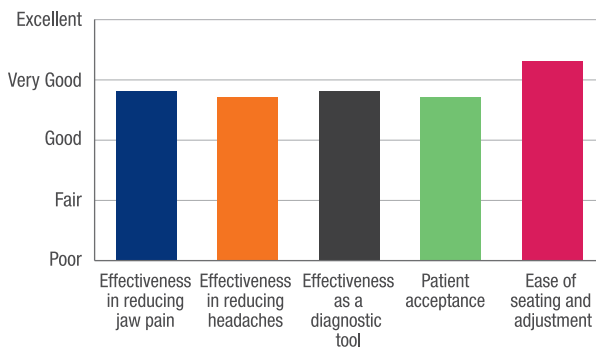

Compared to Competitive Products:

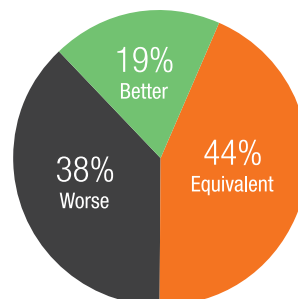

Percentage of Consultants Who Would:

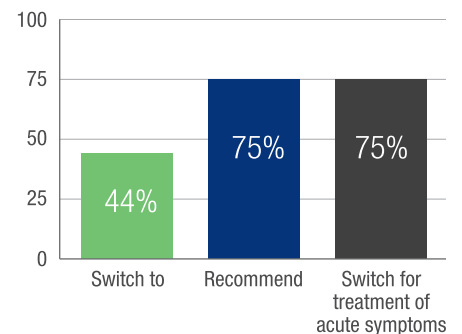

Pulpdent Corporation
www.pulpdent.com

Description

Bioactive restorative materials are a relatively new concept in dentistry. They are resin composites that are reported to release more fluoride than glass ionomers. They also react to pH changes in the mouth and in the process release and uptake calcium, phosphate and fluoride ions, thus helping maintain the chemical integrity of tooth structure. **ACTIVA** is available in two formulations. **ACTIVA BioACTIVE-RESTORATIVE** is more viscous and delivered via a dual-barreled syringe in a delivery gun. **ACTIVA BioACTIVE-BASE/LINER** is less viscous and dispensed via an automix syringe. Both delivery systems are easy to use and designed to be dentist friendly.



Purpose

The purpose of this clinical evaluation was to determine the properties of **ACTIVA BioACTIVE-RESTORATIVE** restorations at placement as well as to monitor its performance clinically at six months.

Clinical Protocol

A total of 117 **ACTIVA BioACTIVE-RESTORATIVE** restorations were placed in 58 patients. These included anterior and posterior restorations (Figure 1) with one, two, three and four surfaces (Figure 2). Sixty-six percent of the restorations were placed with etch followed by the application of a bonding agent. Thirty-four percent of the restorations were placed directly after etching the tooth (without the application of a bonding agent).

The restorations were evaluated at placement for ease of handling, ease of finishing and polishing, ease of adjusting, esthetics, and lack of postoperative sensitivity. Each category was rated on a 1 to 5 scale: 1=poor, 2=fair, 3=good, 4=very good, and 5=excellent.

FIGURE 1 Distribution of **ACTIVA BioACTIVE-RESTORATIVE** Restorations.



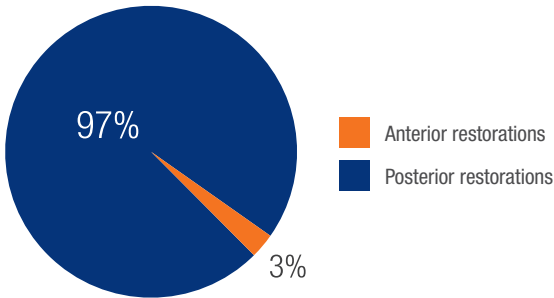


FIGURE 2 Distribution of Surfaces of **ACTIVA BioACTIVE-RESTORATIVE** Restorations.

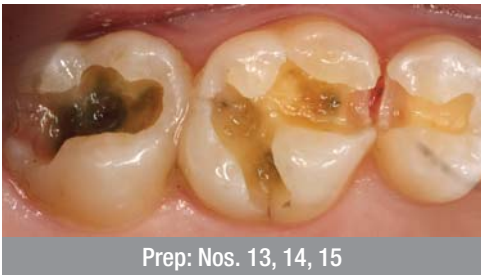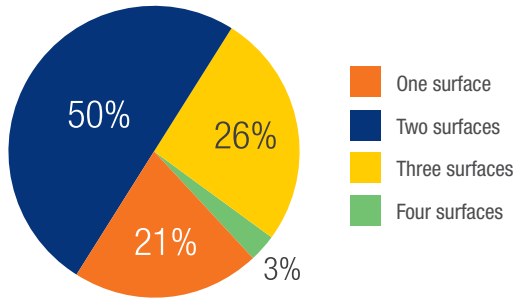

ACTIVA BioACTIVE-RESTORATIVE (cont.)

Clinical Observations at Placement

Ease of Placement

One hundred fifteen **ACTIVA BioACTIVE-RESTORATIVE** restorations received an excellent rating of 5, while two restorations received a rating of 4. In one case, the material was runny. In the second case, the material began to set prematurely.

Ease of Finishing and Polishing

All 117 restorations received an excellent rating of 5 for ease of finishing and polishing.

Ease of Adjusting

One hundred fourteen restorations received an excellent rating of 5 for ease of adjusting while three restorations received a rating of 4.

Esthetics

One hundred fourteen out of 117 restorations received an excellent rating of 5. Three restorations received a rating of 4 because the A3.5 shade looked more like an A2 shade.

Clinical Observations at Six Months

Thirty of the **ACTIVA BioACTIVE-RESTORATIVE** restorations were available for recall at six months. The following parameters were assessed: retention, resistance to marginal discoloration, wear resistance, esthetics, resistance to fracture/chipping, and lack of postoperative sensitivity (as reported by patients). Categories were evaluated on a scale of 1 to 5: 1 = poor, 2 = fair, 3 = good, 4 = very good, 5 = excellent.

Retention

Retention was rated excellent. There were no debonds observed at six months.

Resistance to Marginal Discoloration

Resistance to marginal discoloration was rated excellent. There were no observations of marginal staining at six months.

Wear Resistance

Wear resistance was rated excellent. There were no observations of wear at six months.

Esthetics

Esthetics was rated excellent. The restorations maintained excellent color match with the surrounding tooth structure at six months.

Resistance to fracture/chipping

Resistance to fracture and chipping was rated excellent. There were no observations of fracture or chipping at six months.

Lack of Postoperative Sensitivity

Lack of postoperative sensitivity was rated excellent. Two of the 117 restorations were replaced shortly after placement. In one case a sizable void was noted in the restoration when the x-ray was taken. In the second case the tooth was very sensitive and required root canal treatment. The tooth that required endodontic treatment was a deep restoration and was likely not related to the choice of restorative material. There were no additional reports of postoperative sensitivity at the six-month recall.



No. 14 at 6 months

Summary

Ninety-seven percent of the 117 **ACTIVA BioACTIVE-RESTORATIVE** restorations exhibited excellent esthetics at placement. The ease of placement, ease of finishing and polishing, ease of adjustment, and lack of postoperative sensitivity of the material were rated as excellent. **ACTIVA BioACTIVE-RESTORATIVE** restorations performed well at the six-month recall. THE DENTAL ADVISOR will continue monitoring these restorations over time.

Versa-Temp 2

Sultan Healthcare
www.sultanhealthcare.com




Description

Versa-Temp 2 is a provisional crown and bridge material with increased strength compared to the previous formulation. **Versa-Temp 2** has a natural fluorescence and it has a compressive strength of 450 MPa. It is available in shades BL, A1, A2, and A3. **Versa-Temp 2** reaches a “rubber” phase after 60-70 seconds from initial mix and should be removed from the mouth at that time. It should remain in the matrix until completely cured, 4 minutes from the initiation of mixing. **Versa-Temp 2** comes in 50ml cartridges and fits a standard 1:1 impression gun. **Versa-Temp 2** shade A2 was evaluated by 27 consultants in 406 uses. This provisional material received an 86% clinical rating.

Suggested Retail Cost \$159.99/cartridge

Product Features

Versa-Temp 2 produced strong, durable provisional restorations. It had a relatively fast setting time and was best removed from the mouth during the rubbery phase. No significant shrinkage or distortion occurred, and fit was rated as very good. Margins could be trimmed without chipping, and it polished to a shiny finish. The final esthetics of the A2 shade evaluated were very good, with enough opacity to mask most underlying cements.

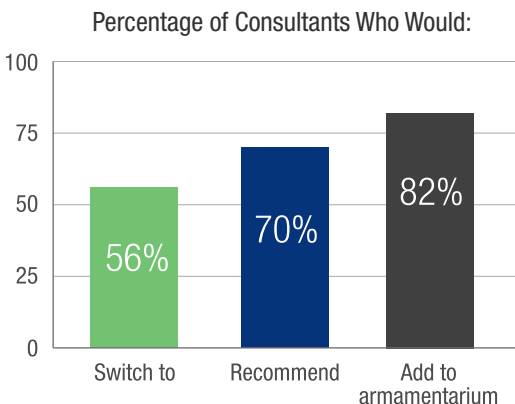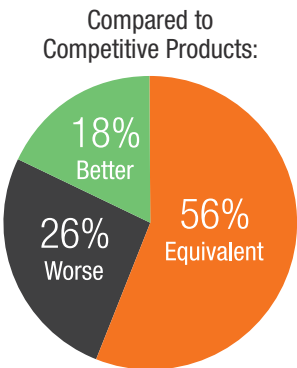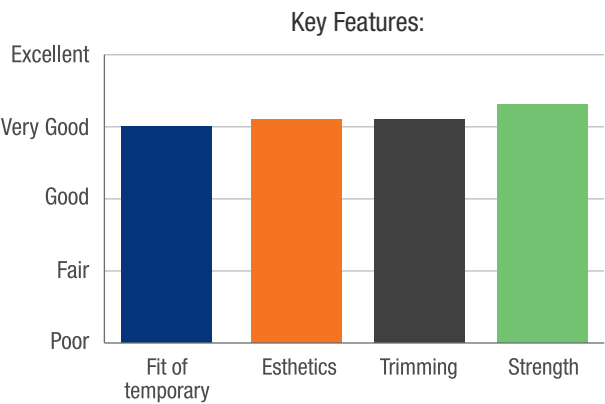

++++



Consultants' Comments

"The strongest temporary material I've used."

"It is was easy to use and held up well during the provisionalization stage of treatment."

"Relatively fast setting time."

"Trims smoothly with a lab carbide."

"Nice esthetics."

Clinical Tip

- Use alcohol to remove the sticky layer after setting. ■

AirBlade AirMirror

UltraLight Optics, Inc
(310) 920-1368 | www.ultralightoptics.com

Description

The **AirBlade AirMirror** is designed to eliminate the need for dental assistants to clean the mouth mirror with an air/water syringe during procedures. The **AirBlade AirMirror** connects to all quick-connect air supply port and uses an air jet to clear water and fog from the surface of the mirror. A thumb knob valve allows the clinician to fully adjust the flow of air. The system comes with four mirrors, four handles and two quick disconnects (quick disconnects are available in all sizes). **AirBlade AirMirror** was evaluated by 13 consultants in 320 uses. This mirror received an 83% clinical rating.

Product Features

AirBlade AirMirror gave more than one third of consultants significant time savings during treatment. Having the mirror constantly clear allows quicker

+++ 1/2



Suggested Retail Cost Complete set of 4 mirror systems \$795

Vera by Young

Young Dental

www.youngdental.com

Description

Vera® by Young™ is a disposable prophylaxis angle designed for improved Visibility, Ergonomics, Reach and Access. The size is small and compact with a slim neck and a smooth, rounded head to help with patient comfort and deliver maneuverability in the mouth in a variety of clinical situations, such as tight buccal mucosa, small mouths and third molars.

Vera by Young is available in two cup styles: the soft, latex-free *Elite Original™* cup with inner and outer ridges for reduced splatter, and improved interproximal cleaning and stain removal; and the soft, latex-free *Petite Web™* cup, which offers flex and flare. **Vera by Young** is available in Classic and Contra styles; **Vera by Young** is available in boxes of 125 or 750 angles. **Vera by Young** was evaluated by 37 consultants in 2095 uses. This disposable prophylaxis angle received a 98% clinical rating.

Suggested Retail Cost

125 angles	\$113.30
750 angles	\$633.45

Product Features

Vera by Young is a solidly built disposable prophylaxis angle that functioned flawlessly during the evaluation period. The fit to the handpiece is excellent, and it runs smoothly with no chatter. The small, rounded

+++++



Consultants' Comments

"I love how I can use every surface of the cup to clean all areas of a tooth."

"Concentric rotation with no wobble."


The slim size is awesome!"

"It is a durable prophylaxis angle that does not break down."


"It is the only prophylaxis angle I have used that it is splatter free."

head and narrow neck provide good access to posterior areas and improved visibility. The ridges on the *Elite Original* cup allow polishing with the inner and outer surfaces of the cup. The ribbing adds stiffness to the cup, and stain removal is very good. ■


Key Features:

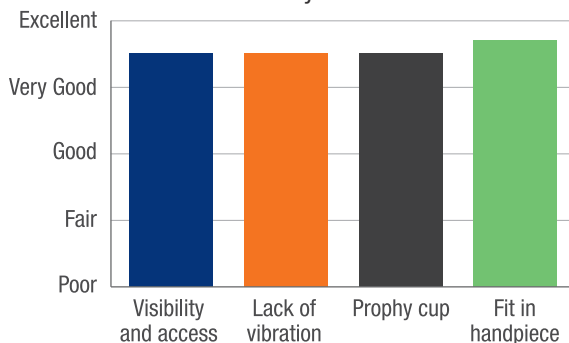

Compared to Competitive Products:

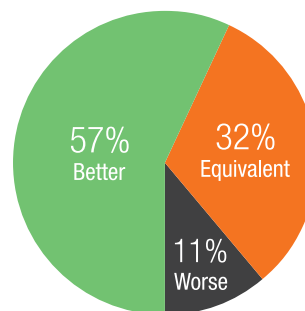


Percentage of Consultants Who Would:

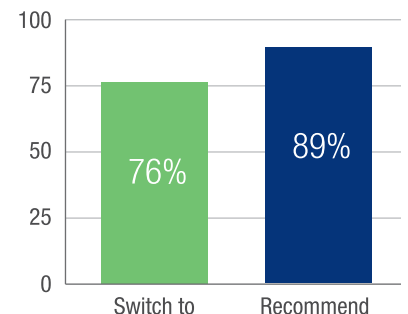

AirBlade AirMirror (cont.)

treatment and frees one of the dental assistant's hands. It was reported that water was cleared effectively from the mirror surface. The airflow knob is sensitive and moves at a light touch to allow full control over the air. Consultants also discovered the airjet cleans the debris and water from the preparation and margins, resulting in extra time savings. The advantage of always having a clear mirror is worth having the lightweight air hose connected to the mirror. ■

Key Features:

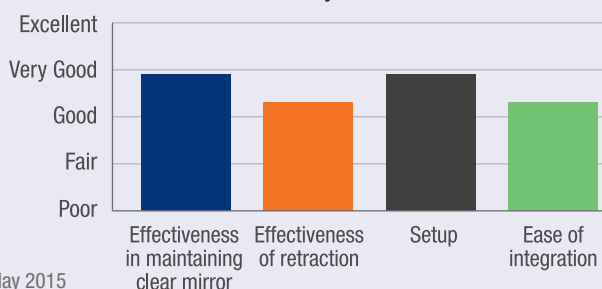


Consultants' Comments

"It does help in keeping the field of vision clear."

"Incredibly helpful - I do not have to stop to wipe the mirror when using the AirBlade AirMirror."

"It is helpful when doing an exam on a mouth breather."

Percentage of Consultants Who Would:

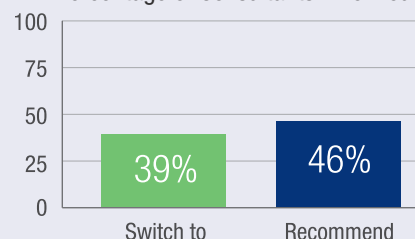

Crystal HD® Mouth Mirrors

Zirc Company

(800) 328-3899 | www.zirc.com

Description

Crystal HD® Mouth Mirrors are designed to help prevent eye strain and stress by displaying clear, shadow-free images. This mirror lens is 40% brighter than traditional rhodium with a 99% true-color reflection. Mirrors come in several styles and sizes: #4 and #5 plastic or stainless heads to fit standard handles and one-piece plastic mirrors with Thin Grip or Soft Grip handles. All mirrors are autoclavable. Plastic handles and heads contain antimicrobial protection to aid in preventing bacterial growth. Plastic heads have a flush fit against the lens, eliminating the gap typically found in metal mirrors and preventing collection of debris. The Crystal HD® #4 Soft Grip Mouth Mirrors evaluated come in a 12 pack of assorted neon colors. The mouth mirrors were evaluated by 33 consultants in 1097 uses. This mouth mirror received a 96% clinical rating.

Suggested Retail Cost

\$74.00 / 12 pack

Product Features

Crystal HD® Mouth Mirrors offer noticeably brighter reflection than standard mirrors. The *Soft Grip* handle is lightweight and ergonomic, but the grooves can trap prophy paste and other debris. The flush design of the mirror head eliminated trapping of debris in this area. Consultants found the plastic material to be more comfortable for patients. Dental assistants rated **Crystal HD® Soft Grip Mouth Mirrors** highly for retraction except in cases of very strong tongues, in which case the plastic handle could bend. The colorful handles were a welcome change from standard stainless instruments.



Consultants' Comments

"The crisp, crystal clear images reduced eye strain!"

"Mirrors were light and pleasant to use."

"Comfortable grip, and patients commented on the bright colors."

"Kid friendly instrument."


"I loved these for viewing the margins on crown preparations."

"Comfort Grip handles are a little shorter than standard handles."


Clinical Tip

- Debris on the mirror face is best removed by wiping with a cleaning solution such as the Mirror Magic System. ■


Key Features:

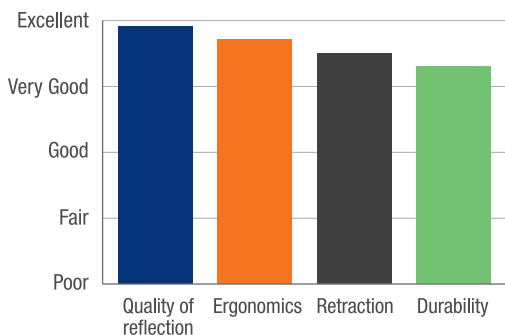

Compared to Competitive Products:

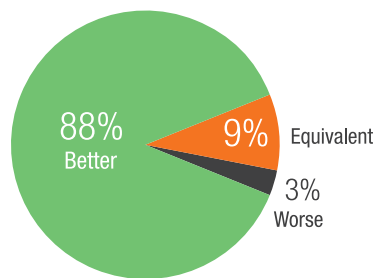

Percentage of Consultants Who Would:

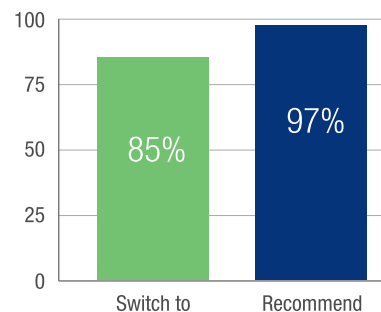

SPECIAL THANKS TO:

Select Senior Clinical Consultants (Over 20 years):

J. Amara, CT · R. Fisher, OH · W. Gregory, MI · E. Katkow, MD · J. Lockwood, MI · J. Mayer, OH · W. Nagy, TX · G. Poy, MI · J. Shamraj, MI · R. Trushkowsky, NY · P. Yaman, MI

Senior Clinical Consultants (15-19 years):

S. Harlock, MI · R. Anthony, OH · K. Baker, TX · F. Berman, PA · L. Brimhall, MT · M. Briskin, NY · L. Brown, MI · R. Ciccone, MI · C. Colbert, MI · M. Conrad, PA · J. Dingman, WI · J. Doueck, NY · M. Dwoskin, MI · M. Eannaccone, NY · R. Engle, IL · K. Fairbanks, MI · K. Fischer, IN · G. Franco, NY · N. Garlisi, OH · K. Goodman, MI · S. Graber, IL · P. Grandsire, NY · B. Gursky, MI · E. Gutman, NY · K. Hamlett, TX · G. Hart, OH · R. Herwig, KS · J. Kaminski, MI · R. Kaprielian, NJ · M. Kastner, OH · C. Kehr, MI · M. LaMarche, WA · J. Leitner, MI · R. Lezell, MI · M. Livernois, NC · B. Manne, FL · N. Mansour, MI · N. Markarian, CA · J.W. Mikesell, IL · G. Mosso, PA · E. Mosso, PA · E. Odenweller, OH · J. Paris, TX · D. Parris, GA · D. Peterson, MD · T. Pieper, WY · D. Pitak, MI · V. Plaisted, NY · D. Qualliotine, NC · C. Reed, MI · P. Sandvick, WI · K. Schwartz, FL · J. Shea, MO · B. Shumaker, NJ · B. Sims, NY · H. Tetelman, OH · C. Trubschenck, CA · P. Tu, CA · S. Ura, NH · W. Walcott, MI · L. Wee, MI · D. Wojtowicz, MI · M. Zahn, MI

Clinical Consultants (14 years or less):

A. Albright, NY · N. Alexa, MI · G. Alfe, IL · B. Argersinger, NC · P. Arsenault, MA · A. Bagchi, MI · B. Barricklow, OH · L. Bartoszewicz, MI · J. Bechtel, MI · C. Bhatti, MI · L. Bishop, MI · T. Bizga, OH · G. Bloomfield, MI · J. Bostic, OH · C. Brown, IL · W. Brownscombe, MI · E. Brust, MI · S. Bunek, MI · J. Bunek, MI · T. Burns, MI · J. Bush, PA · H. Cadorette, MI · P. Campo, NY · P. Cracchiolo, MI · D. Chacko, TN · P. Chaiken, IL · D. Chenevert, MI · W. Co, MI · M. Connelly, MI · C. Dietz, OH · S. Dillingham, NY · R. Dost, VA · A. Dutko, MI · M. Eiford, MI · O. Erdt, MI · F. Facchini, MI · F. Falcao, FL · M. Feinberg, NY · L. Feldman, NJ · M. Frankman, SD · M. Glovis, MI · C. Goldin, MI · M. Grant, MI · B. Greenwood, UT · J. Griffin Jr., MO · P. Gronet, MI · H. Gulati, MA · D. Haas, Ontario · J. Haddad, MI · S. Harden, GA · J. Hastings, CA · D. Hauck, CA · J. Hengehold, MI · B. Herrmann, CA · D. Hock, MI · C. Huang, CA · P. Indianer, MI · J. Ireland, MI · C. Jaghab, MI · J. Jaghab, MI · M. Kachi-George, MI · D. Kapp, NY · J. Karam, MI · E. Kelly, GA · J. Kelly, GA · L. Kemmet, MN · D. Keren, NY · C. Knapp, MI · M. Koczarski, WA · B. Kolb, MI · D. Kuras, MI · C. Latham, OH · R. Le, NC · S. Lever, MD · I. Levine, NY · E. Lowe, BC, CAN · J. Lusby, MI · M. Man, NY · C. Manduzzi, MI · K. Mantzikos, NY · F. Margolis, IL · M. Mason, NY · T. McDonald, GA · C. McLaren, MI · J. McLaren, MI · M. McMullin, MI · H. Menchel, FL · M. Migdal, MI · J. Minsky, CA · R. Mizrahi, NY · M. Moeller, MI · G. Molinari, MI · L. Montes, NY · M. Murphy, MI · J. Nash, MI · A. Nazarian, MI · N. Nealis, IL · J. Olitsky, FL · J. Olsen, MI · F. Orlando, NY · R. Oshrain, NY · A. Overall, MI · S. Owens, MI · M. Paquette, MI · J. Park, IL · J. Parrott, MI · M. Patel, MI · C. Pike, MI · B. Pittsley, MI · T. Poirier, MI · D. Radtke, MI · G. Raichelson, Ontario · G. Ramos, NY · C. Ramsey, FL · T. Reeves, TN · N. Rego, CA · G. Reskakakis, NY · J. Riggs, MI · J. Rowe, AR · A. Saddy, MI · L. Seluk, MI · R. Selvan, NJ · M. Shapiro, MI · A. Shemesh, IN · P. Shumaker, MI · B. Silver, MI · S. Simpson, MI · J. Smith, MI · C. Stevens, OK · B. Stevenson, MI · B. Stieper, MI · S. Tamber, MI · G. Tarantola, FL · T. Teel, IN · Lyn. Vandelaar, MI · H. Vann, MS · M. Waranowicz, MI · H. Whitt, MI · L. Williams, MI · K. Wilson, MI · D. Wolf, MA · W. Wright, CA · N. Yeager, MI · H. Yeung, CA · D. Young, MI · P. Zanetti, MI · J. Zanetti, MI · S. Zimmer, MI

Laboratory Consultants:

Apex Dental Milling, MI · Bullinger Dental Lab, MI · Centric Dental Laboratory, TX · Cornerstone Dental Studio, Inc., MI · David's Dental Laboratory, NY · Expertec Dental Lab, MI · Heritage Dental Laboratory, IL · Technique Crown and Bridge, Inc., NC